

SPA safety Screwdriver Kit

Reduce the Risk of
Accidental Swallowing



Instagram



Facebook

Comfort with Moderation
Invisible Protection




Code	Compatibility
HEX-1.00	Ankylos
HEX-1.16	Anthogyr
HEX-1.20	Osstem、Biomet 3i、DIO-UF、Megagen、NEO-Biotech CWM/INNO、Evertis、YES、Anker、T-Plus、SIC、SG Datsing、xive、Leadoss、CSM、Zenois、S.I.N ICX/ACM-3.3、WARANTEC、Biotech Dental OneDayBioTech、DentalMax-DM、Trausim-RS
HEX-1.25	MIS/Lance-seven、BEGO、Alpha-BioTech、Biodenta、AB ANGELS、Adin-Touareg、Implant Direct Biohorizons Cortex、Dentin、DENTIS、Martina、DentalMax-DD、IEZ
HEX-1.27	Dentium-superline/simperline、Ticare、BAM、Medentika MIS-C1、B&B/C-tech、Camlog/Conelog、Trausim-RD、CLC
HEX-1.30	Astra、WEGO
S-TORX	Schroeder、CANSUN/TAISHAN、Trausim-BL/BLE ITI Straumann、CRS、Bioconcept-BC
N-TORX	Nobel BioCare、Trausim-NT、Dental Master-DM

SPA Safety Screwdriver Kit

Accidental Drop, Instant Protection | Reduce Surgical Risks

Different Colors for Identification Easy to Organize

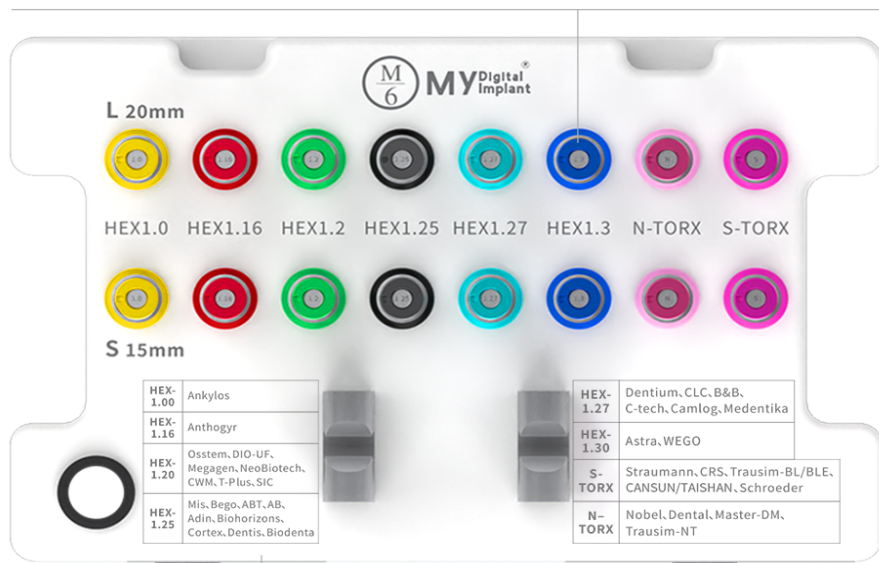

8 Head Types 2 Lengths Multiple Sizes

Designed for 99% Clinical Cases

SPA Safety Screwdriver

Reduce the Risk of Accidental Swallowing
Medical-Grade Stainles



Top Laser Marking


Quick Size Identification



Strong Retention

Perfect-Match with
Screw Interface

Secure Locking, No Loosening

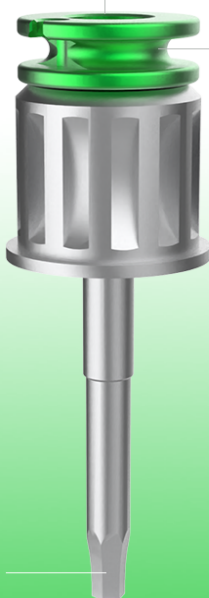

**Safety String
Retention groove**

Designed for Safety
String to Fix Securely

SPA Safety String

Simple Application | Convenient to Operate | Cost-Effective

A Simple Solution for Drop Prevention

Elastic Retention

Easy to Use

100PCS/Box

Safety String
Prevention Approach

Silicone
Elastic Retention

Easy-Open
Individual Packaging

100 pcs/box

SPA Safety Screwdriver and String Kit Contents

SPA Safety String(100PCS/Box)



SPA Safety Screwdriver(16 Screwdrivers)

L
20mm



S
15mm



SPA Safety Screwdriver and String How to Use

Select



1. Open the tool kit and select the appropriate screwdriver based on the laser marking or color

Place



2. Secure the safety string to the retention groove on the screwdriver

Use





3. Place the other end of the safety string around your finger

Reduce the Risk of Screwdriver Dropping

Help Ensure a Successful Procedure

Accidental Dropping of the Screwdriver into the Airway
Real Clinical Case Extremely High Risk

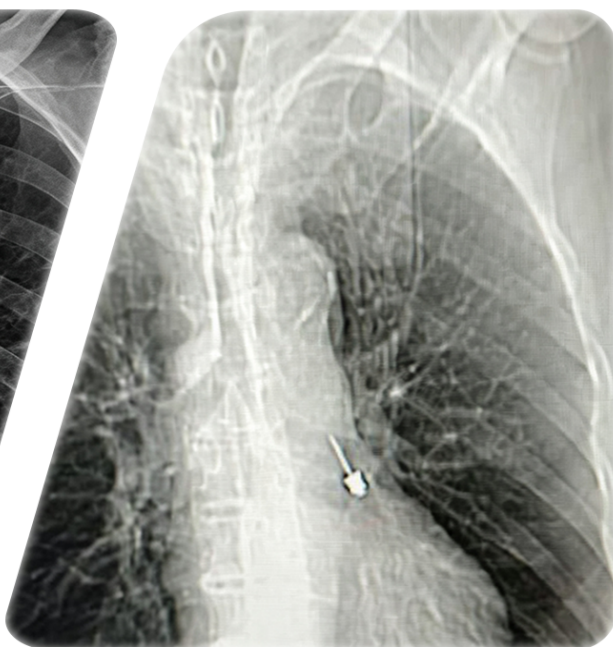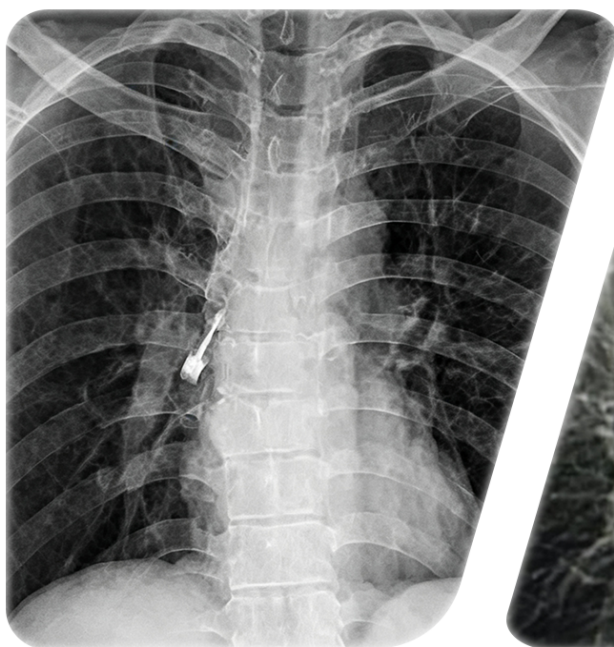

Sourced: Yangzi News